


Conversor TV Pro

The Clarity of Sound

Assistive Listening Devices

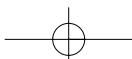

Easy Listening

Introducing The Conversor® TV Pro – an affordable assistive listening device that offers users the versatility to enjoy the most out of TV and audio listening.

Whether at home watching the TV, listening to the radio, or in a classroom, conference or church hall, lecture or meeting room, hearing impaired participants face a real challenge understanding the speaker clearly amongst surrounding background noises.

Conversor TV Pro – is for the individual user to listen to the TV or other Audio source without increasing the volume on the TV which can irritate other viewers. Simply plug the transmitter into the TV, place the Receiver around the neck, switch the hearing aid to Telecoil and listen through the hearing aid or just plug in a set of headphones.



Conversor TV Pro – addresses the needs of a hearing impaired participant, whether a Hearing Instrument wearer or for use with headphones or earbuds, in classroom, conference or church hall, lecture and meeting room settings, providing front row sound level and clarity up to 50 metres away from the speaker. Just plug the TV Pro into the Sound System and it broadcasts clearly to the Receiver worn around the participant's neck.

Conversor TV Pro – includes one Conversor TV Pro Transmitter which is connected to the Audio or TV output of the Sound System and a Pendant Receiver. The Transmitter can be used with multiple Receivers which are set to the same frequency. The TV Pro Transmitter is mains powered and the Receiver can be charged using the same charger.

A Wireless link to the hearing Instrument – the Conversor TV Pro Transmitter, which has an automatic gain control to maximise sound quality, broadcasts the sound directly to both hearing instrument wearers and headphone or earphone wearers enabling them to hear the television or sound source clearly. The desired sound is transmitted with enhanced clarity to the Receiver (worn around the user's neck either above or below clothing), which sends the signal to the hearing instrument via the inductive neck loop or

directly to the headphones or earphones. By transmitting directly there is no interference from background noises and the desired sound is much clearer and intelligible, even in noisy classroom or meeting situations. By increasing the available sound level, less strain is placed on hearing instrument amplifiers, eliminating distortion and feedback.

TV/Conference Pro Multipack – is totally portable and enables Schools, universities, churches and conference venues to fully comply with the DDA requirements for the Hearing Impaired. It avoids the expense of installing a fixed room loop and for groups that require mobility this can be the ideal complement to a Sound System.

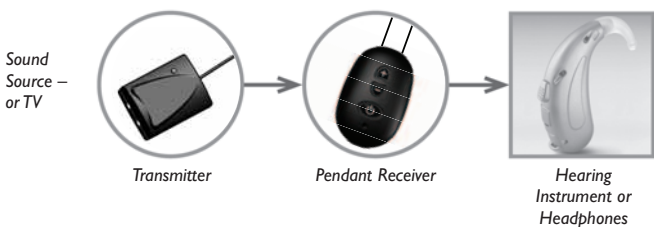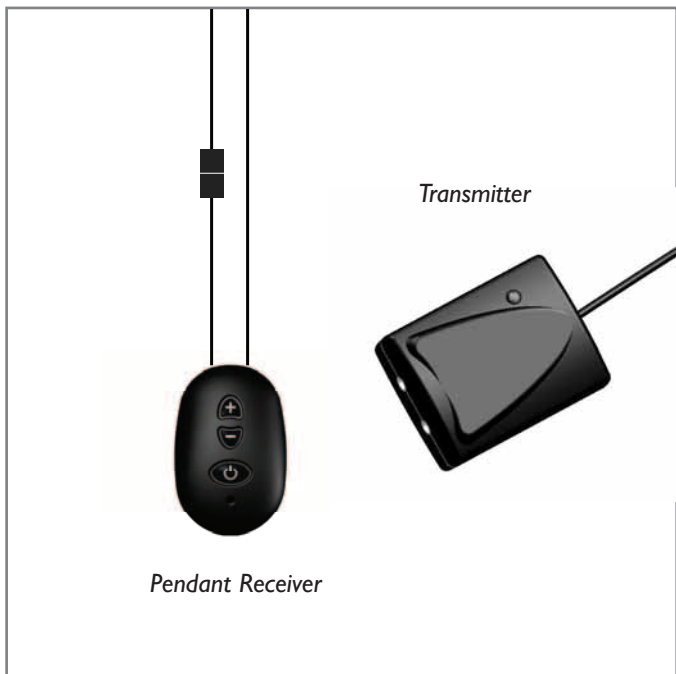

Technical Specifications:

- Input: 100 – 240 Vac
- Output: 5Vdc 2.6 Amp
- Charge Time: 2.5 Hours
- Battery Life: up to 8 Hours
- Audio input levels accepted between 50mV and 5V RMS automatically set, so Pro may be connected to Line, headphone or speaker-level signals.
- Phono, SCART, 2.5mm or 3.5mm stereo jack connector adapters provided for audio input



Conversor Limited

The Lansbury Estate,
102, Lower Guildford Road,
Knaphill, Woking, Surrey,
GU21 2EP (UK)

T: +44 (0) 870 066 3499

Freephone: 0800 298 9982 (UK only)

F: +44 (0) 870 066 3669

E: info@conversorproducts.com

W: www.conversorproducts.com

TV Pro BRO